

Canine Rehab Systems


Specifications Sheet

Portable Cavaletti System

Model PCS-03

GENERAL:

Our Model PCS-03 Portable Cavaletti System is designed for the mobile practitioner or small clinic, where space and portability are important. The frame folds down to a size which will fit in the trunk or back seat of almost any vehicle, and is able to be set up in minutes for use. The rails are individually adjustable in height, allowing for use with all sized dogs, and can even be off-set side to side, for additional proprioceptive challenges.

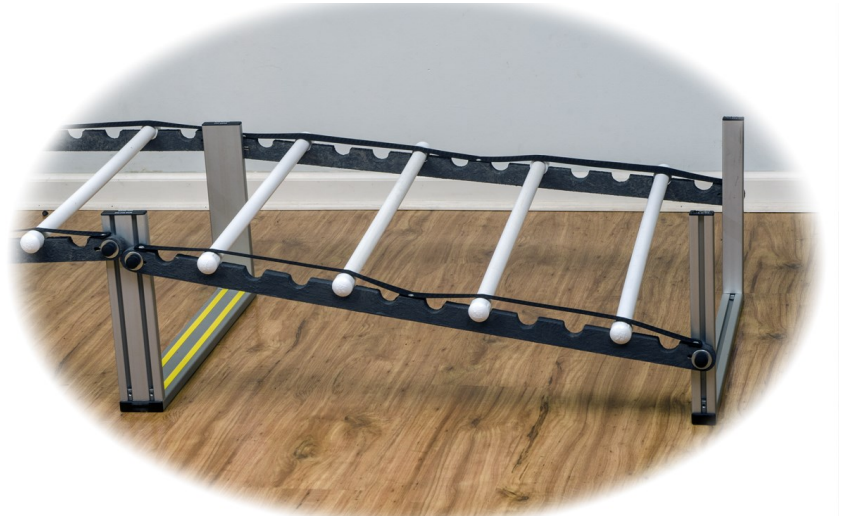

DESIGN FEATURES:

Model PCS-03 is constructed using 3/4 in aluminum extrusion uprights for strength, durability, and dimensional stability. Side rails are machined hardwood. Overall working length is 84 in, with a width of 26 in. Folded for storage or transport, the PCS-03 has a length of 43 in.

Each rail is independently adjustable by use of an end mounted locking knob. Rail height may be set between 2.5 in and 12 in. Rail supports are spaced at 3 in intervals. Heavy duty woven elastic is attached at intervals to secure rails in place.

A set of eight, 3/4in PVC rails with capped ends are included, with a cloth bag for storage / transport.

Canine Rehab Systems LLC

CanineRehabSystems.Com



3787 Jacks Creek Rd NW
Monroe, GA 30655
Info@CanineRehabSystems.Com